


PRODUCT INFORMATION

# MODUS 2 Midface



**MODUS**

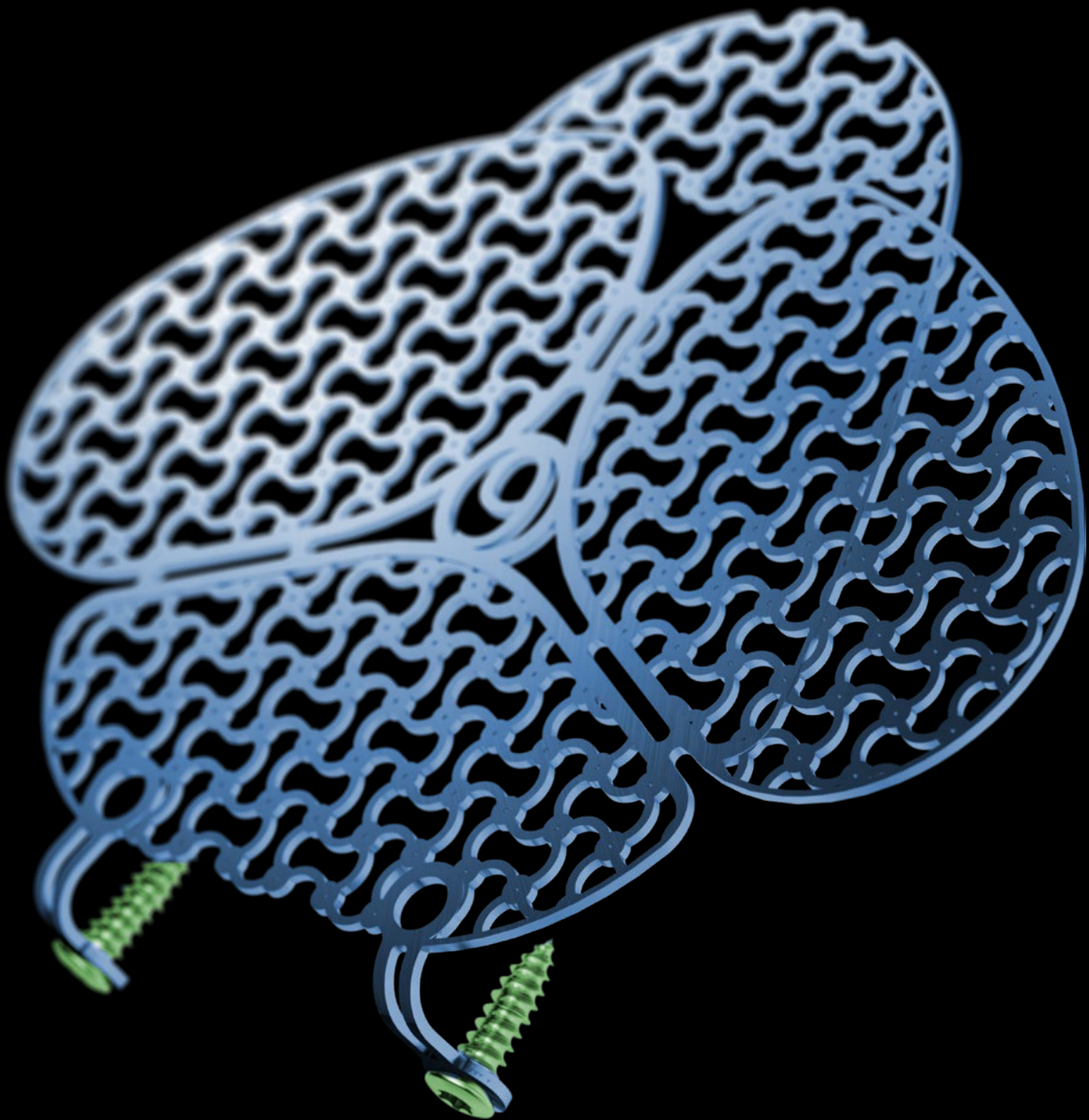



# Contents


- 4 Explore Various Perspectives
- 5 Plates
- 6 Orbital Plating System OPS
- 7 Mesh Plates
- 8 Technology
- 9 Cortical Screw Design
- 10 Instruments
  - Orbital Retractors
  - Mesh Cutting Pliers
- 11 Storage
- 12 Ordering Information

For further information regarding the MODUS 2 product line visit [www.medartis.com](http://www.medartis.com).

# Explore Various Perspectives

The MODUS 2 Midface system is designed to support surgical workflows with a wide selection of plate designs and a unified instrumentation approach. All screws are compatible across plates, and a single screwdriver is used throughout the system, reducing the number of required instruments.

An intuitive, color-coded identification system for twist drills and screws facilitates selection during procedures. The system also includes a comprehensive screw portfolio, such as self-drilling SpeedTip screws, and is available in sterile packaging to support procedural organization and efficient handling.

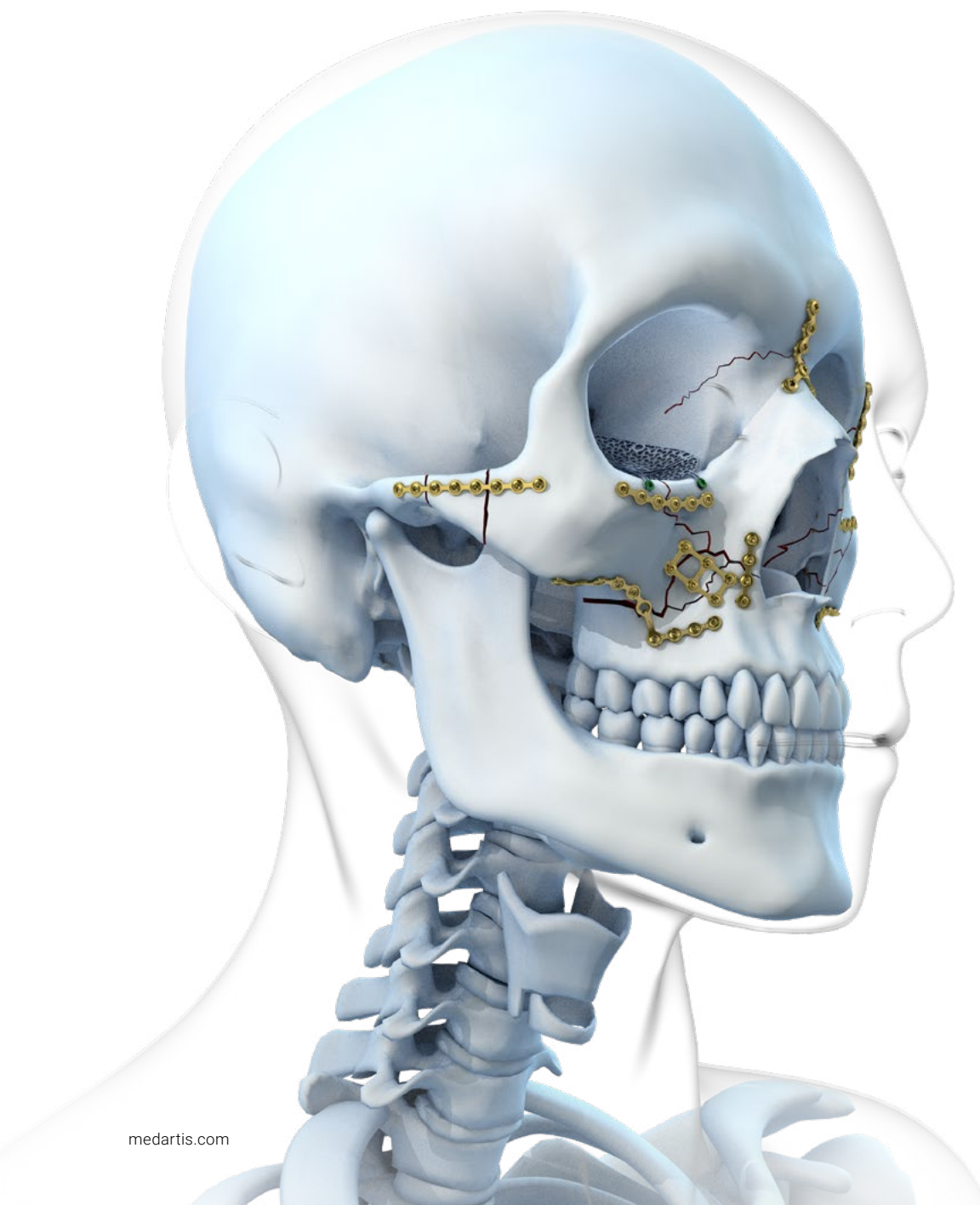


# Plates

MODUS 2 Midface plates offer a broad selection of anatomical and indication-specific plate designs to support midface trauma management, fixation of maxillary osteotomies, and reconstructive procedures. The system is organized by plate thickness, ranging from 0.5 mm to

1 mm, offering versatile options tailored to diverse clinical indications and surgeon preferences. Developed to support intraoperative flexibility and precise adaptation, MODUS 2 Midface plates are a practical solution for midface procedures.

A range of plate thicknesses to accommodate all anatomical midface regions in terms of bone thickness, soft tissue coverage and stability

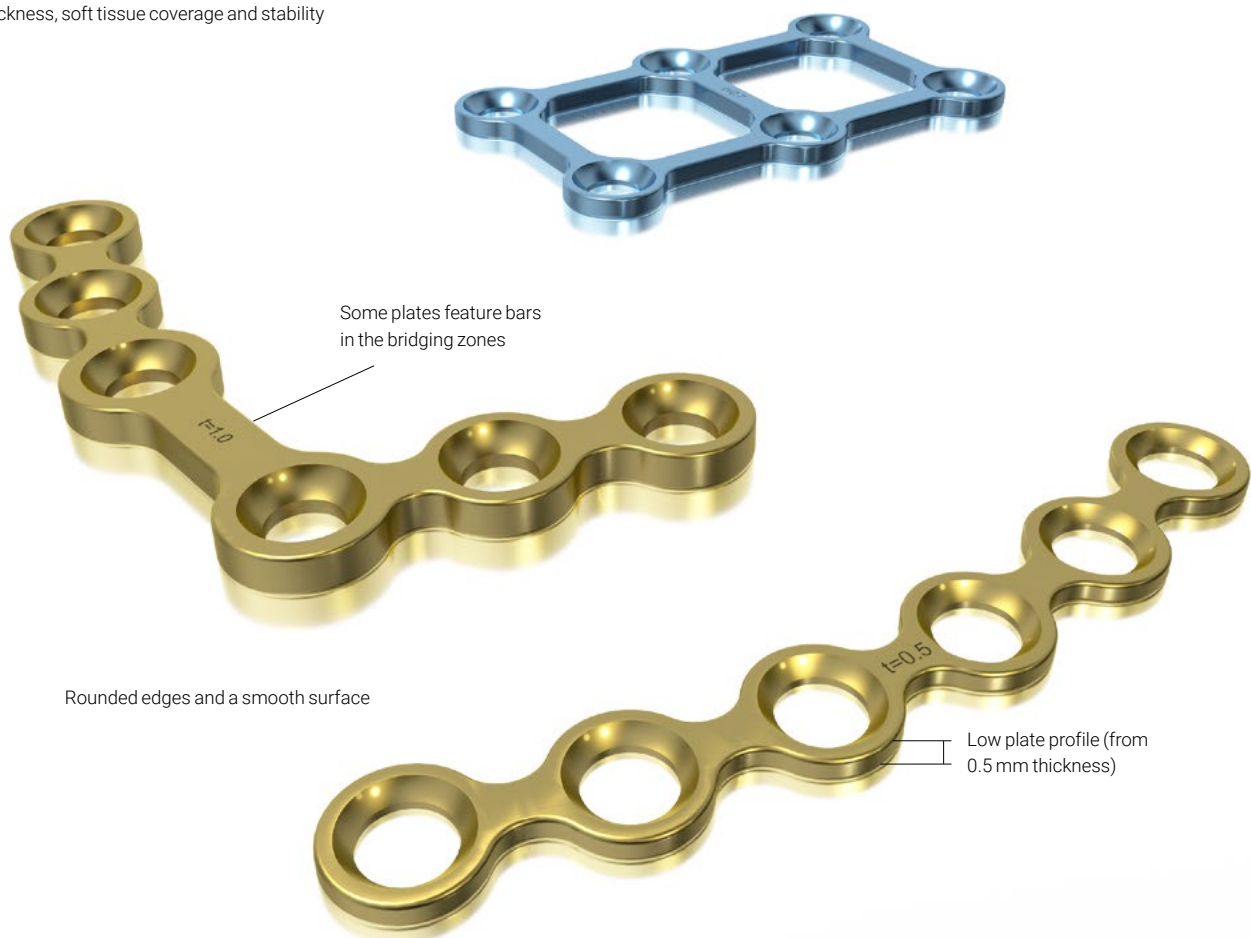

# Orbital Plating System OPS


## Adaptable solution for orbital reconstruction.

MODUS 2 Orbital Plates are designed to support precise reconstruction of orbital floor and wall defects.\*

Manufactured from flexible pure titanium, the implants feature low profile heights of 0.35 mm and 0.4 mm to support precise contouring and anatomical adaptation. The

3D mesh structure reflects orbital anatomy and includes predefined cutting zones. Integrated bars add structural stability where needed.

Available in four implant geometries, the system is designed to accommodate a range of fracture patterns across categories I–IV.

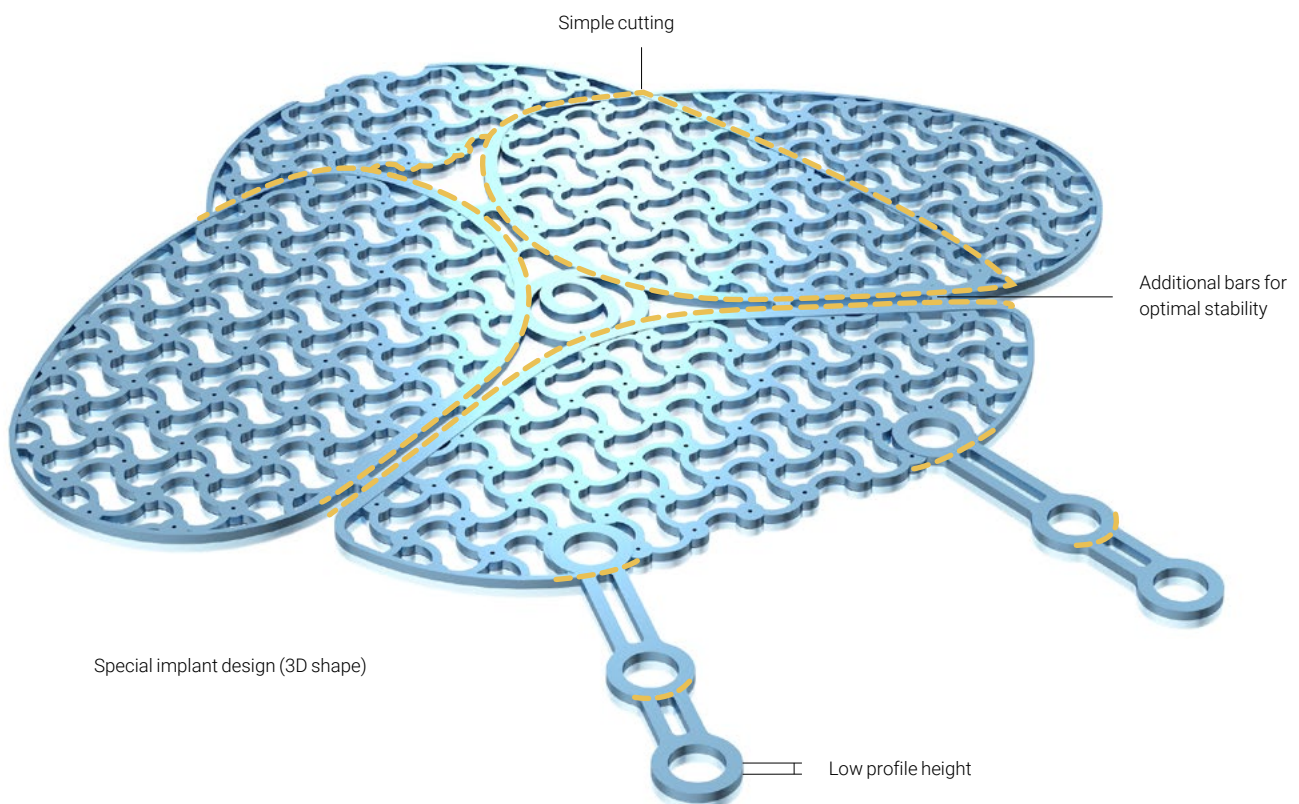

\* Jaquiéry C, Aeppli C, Cornelius P, Palmowsky A, Kunz C, Hammer B. Reconstruction of orbital wall defects: critical review of 72 patients. Int J Oral Maxillofac Surg. 2007;36(3):193-199. doi:10.1016/j.ijom.2006.11.002

# Mesh Plates

## Versatile design.

MODUS 2 Mesh Plates are available in two sizes, each offered in three different thicknesses, 0.25 mm, 0.4 mm, and 0.6 mm, to accommodate a broad range of clinical needs. The 3D mesh structure is engineered for adaptability, with a design that allows for screw head countersinking when needed.

Each plate features two adjacent rounded corners (50 × 50 mm size) to facilitate positioning, along with smooth surfaces and rounded edges.


3D mesh design

The plate design allows countersinking of the screw head

Low profile height:  
0.25 mm / 0.4 mm / 0.6 mm

Rounded edges and a smooth surface


Two adjacent rounded corners

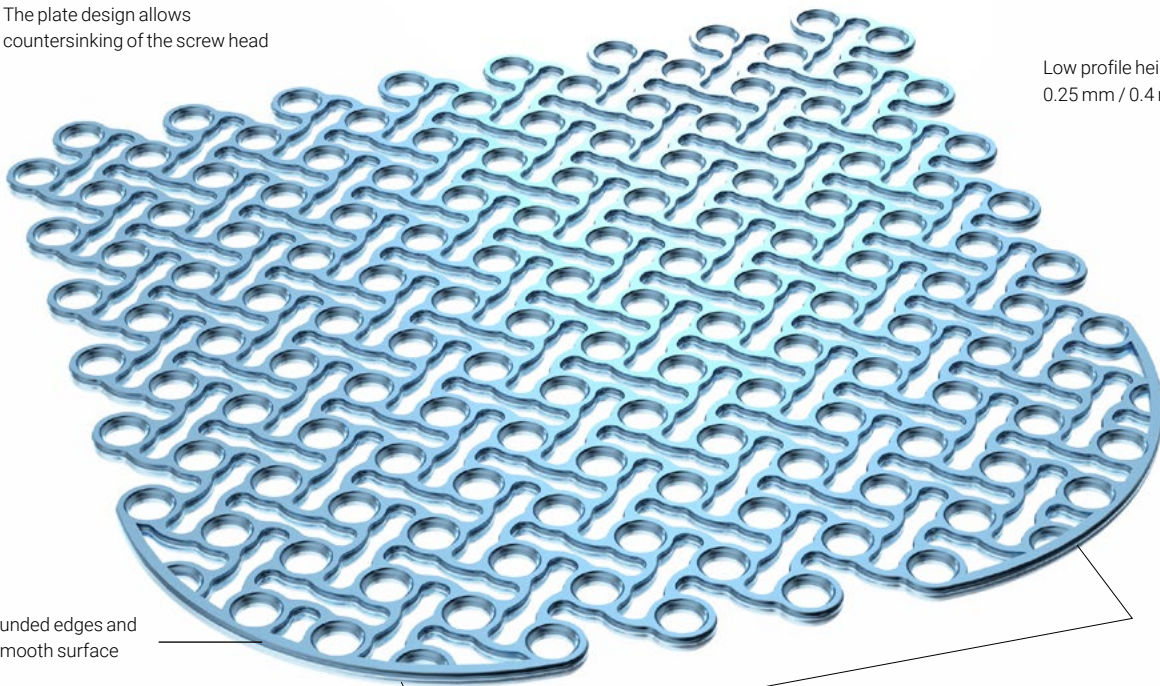

# Technology


## Screw Features

- HexaDrive screw head design:
  - HexaDrive interface with self-holding properties between screw and screwdriver
  - Increased torque transmission
  - Simplified screw pick-up due to the self-holding technology
- Atraumatic screw tip provides soft tissue protection when inserting screws bicortically
- Soft tissue protection due to smooth screw head design
- Double-threaded screws reduce screw insertion time
- Increased torsional, bending and shear stability due to conical core
- Precision cut thread profile for sharpness and self-tapping properties

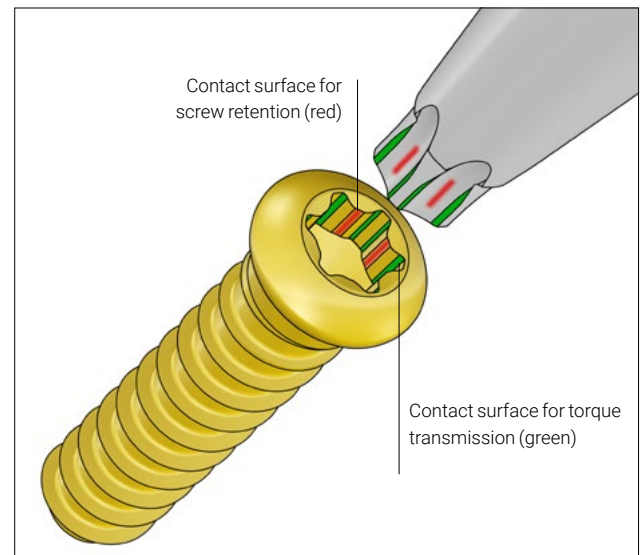

## SpeedTip Thread Design



- Functionally unique cutting with immediate bite
- Immediate cutting of the bone with only slight axial pressure
- The triangular tip design permits simultaneous drilling, tapping, and compression of the bone tissue\*
- The triangular tip design results in increased pull-stability\*
- Reduced insertion torque thanks to the polygonal tip and tapered shaft

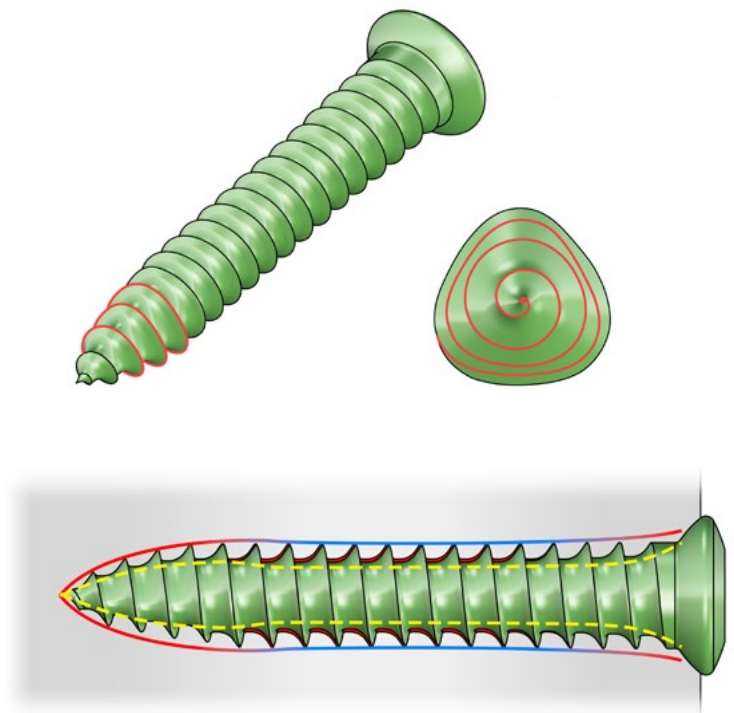

\* Heidemann W, Terheyden H, Gerlach KL. Analysis of the osseous/metal interface of drill free screws and self-tapping screws. J Craniomaxillofac Surg. 2001;29(2):69-74.

# Cortical Screw Design


MODUS 2 cortical screws feature an optimized thread design at the screw head and a tapered shaft, enhancing insertion behavior.

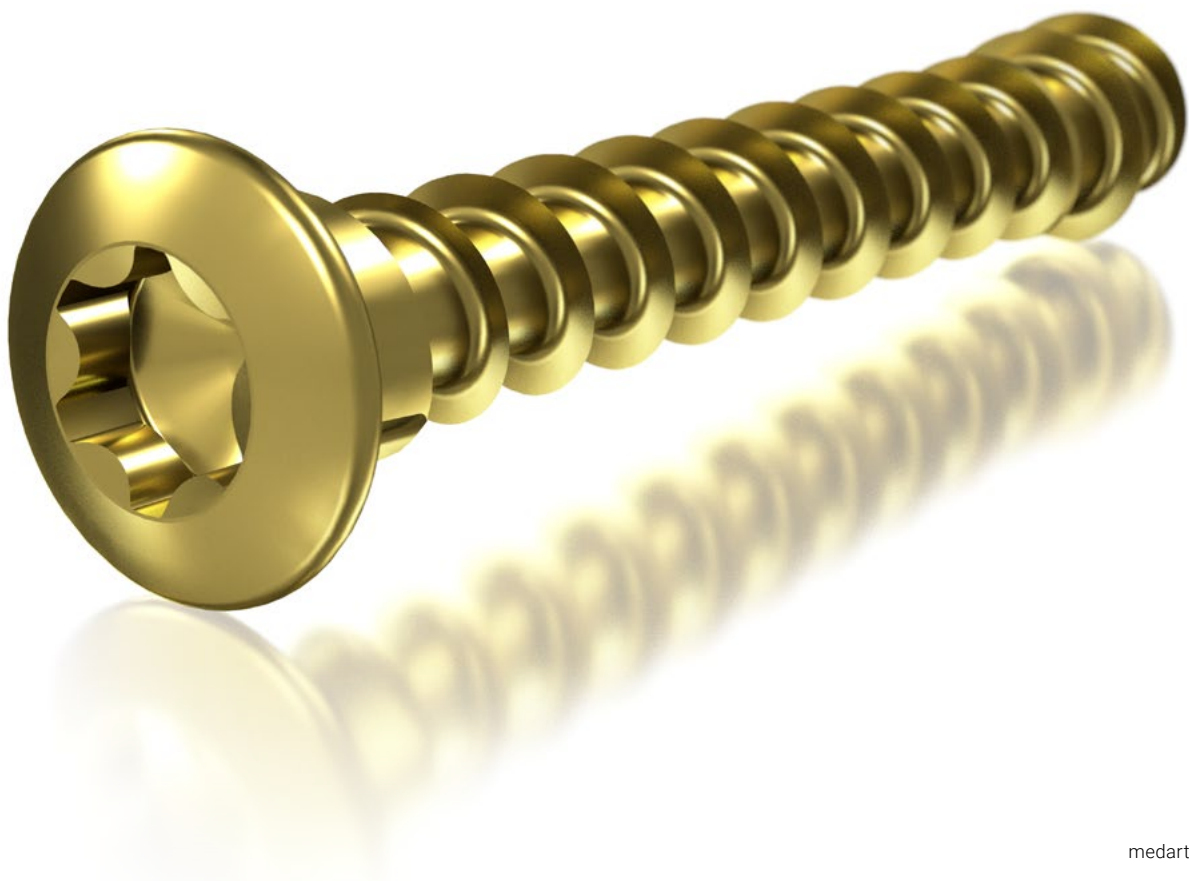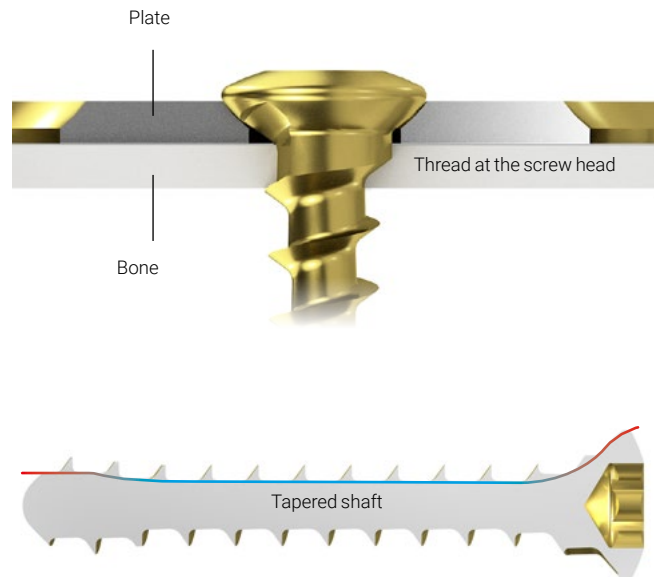

# Instruments

## Orbital Retractors

The MODUS 2 Midface OPS includes three orbital retractors – left, right, and neutral – to protect orbital soft tissue and to determine the defect size. Each retractor features both small and large ends, with shapes and dimensions specifically matched to the MODUS 2 Midface orbital floor plate designs for optimal fit and access.

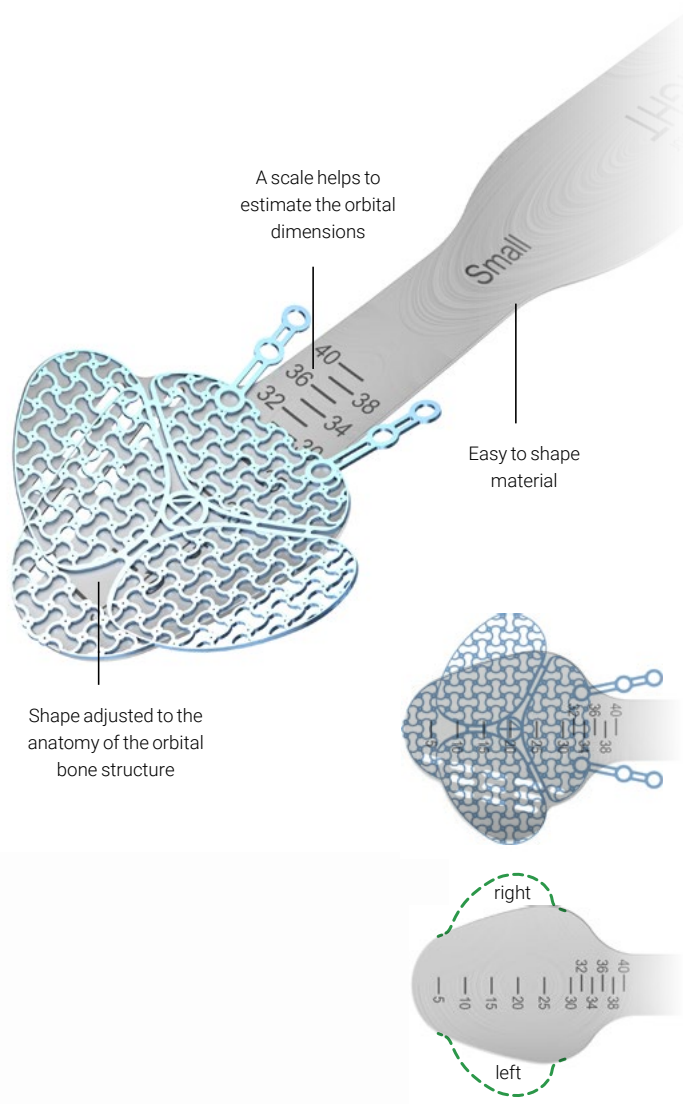


## Mesh Cutting Pliers

The mesh cutting pliers can cut any type of mesh. The cutting process rounds off the cut edge.



# Storage

- Compact, clearly arranged system
- Easy to use
- User-friendly storage of implants and instruments
- Validated cleaning and sterilization of the implants in the container

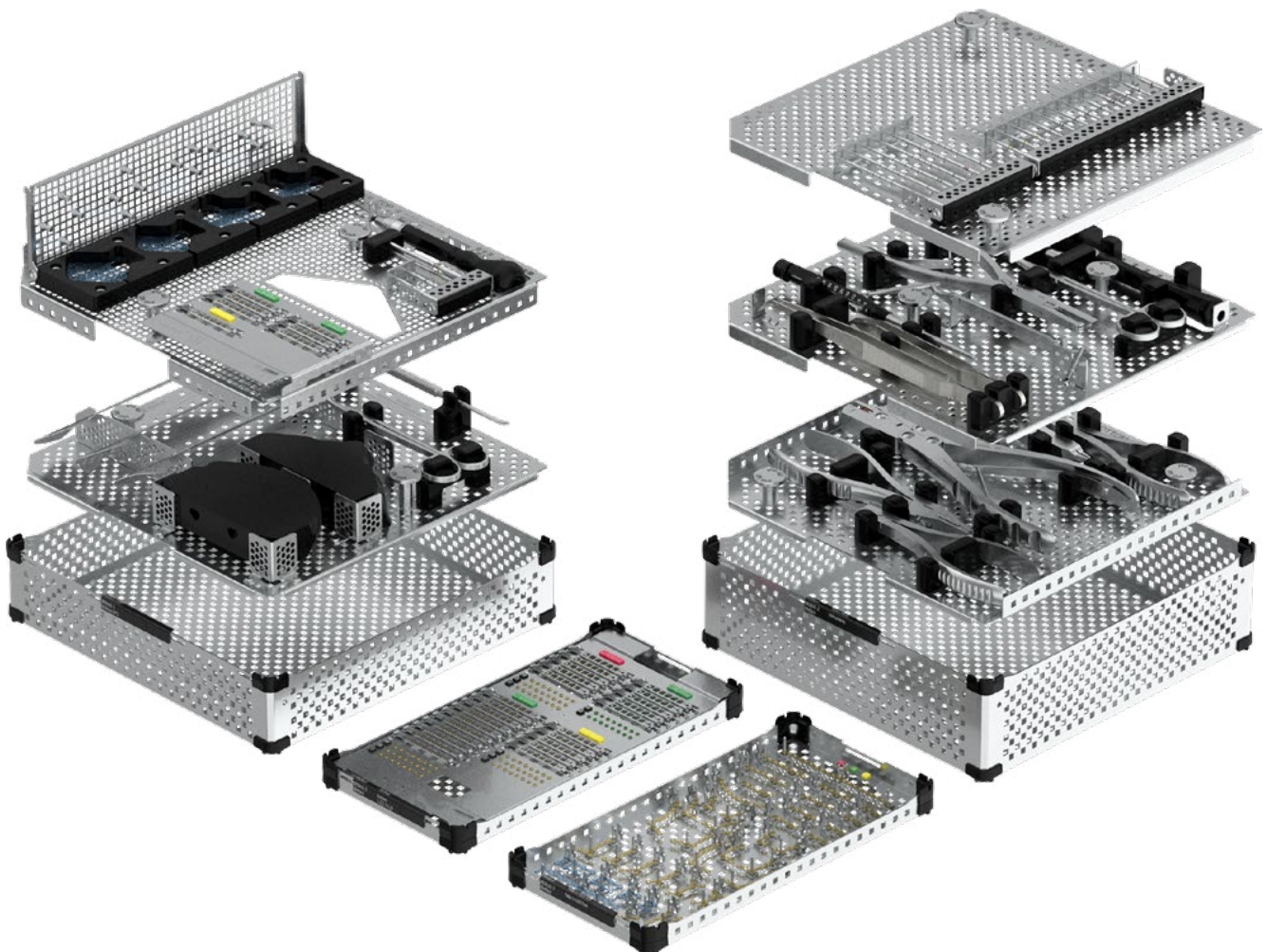


# Ordering Information

## Orbital Floor Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67), semi-rigid  
Plate thickness: 0.35 mm / 0.4 mm

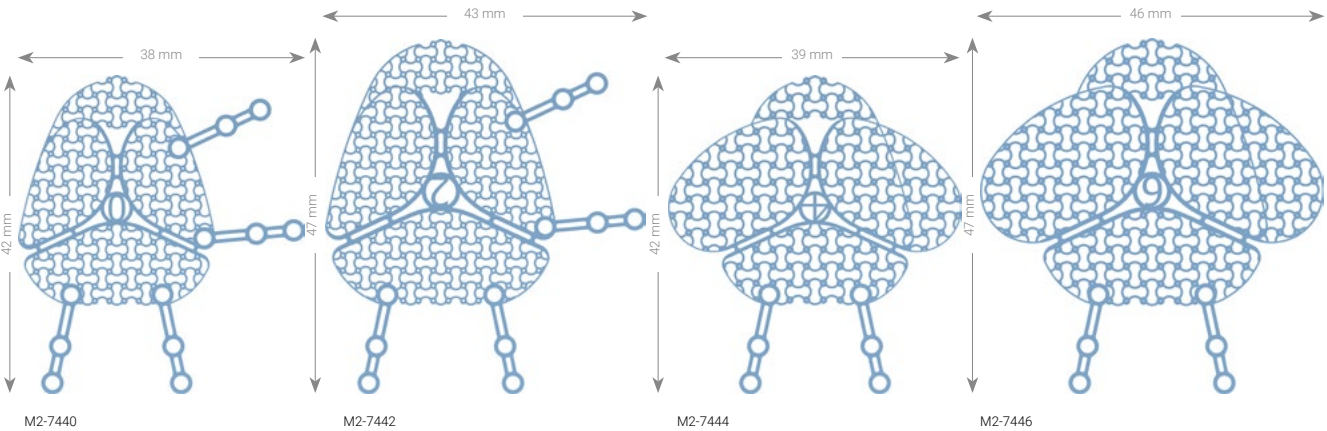


| Art. No. | STERILE  | Description | Plate Thickness | Holes | Pieces / Pkg |
|----------|----------|-------------|-----------------|-------|--------------|
| M2-7440  | M2-7440S | 1 – 2 small | 0.35 mm         | 12    | 1            |
| M2-7442  | M2-7442S | 1 – 2 large | 0.35 mm         | 12    | 1            |
| M2-7444  | M2-7444S | 2 – 4 small | 0.4 mm          | 6     | 1            |
| M2-7446  | M2-7446S | 2 – 4 large | 0.4 mm          | 6     | 1            |

Scale 1:1

## Mesh

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67), semi-rigid  
Plate thickness: 0.25 mm / 0.4 mm / 0.6 mm

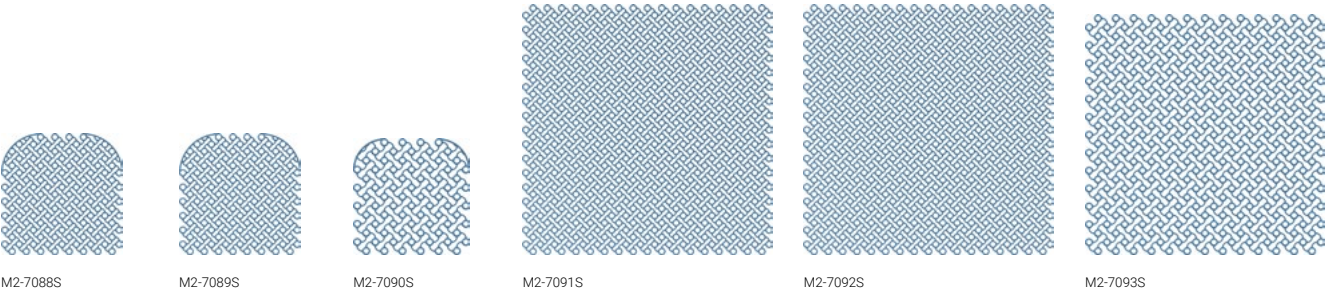


| Art. No. | STERILE  | Description | Size         | Plate Thickness | Pieces / Pkg |
|----------|----------|-------------|--------------|-----------------|--------------|
| M2-7088S | M2-7088S | mesh 1      | 50 × 50 mm   | 0.25 mm         | 1            |
| M2-7089S | M2-7089S | mesh 2      | 50 × 50 mm   | 0.4 mm          | 1            |
| M2-7090S | M2-7090S | mesh 3      | 50 × 50 mm   | 0.6 mm          | 1            |
| M2-7091S | M2-7091S | mesh 1      | 100 × 100 mm | 0.25 mm         | 1            |
| M2-7092S | M2-7092S | mesh 2      | 100 × 100 mm | 0.4 mm          | 1            |
| M2-7093S | M2-7093S | mesh 3      | 100 × 100 mm | 0.6 mm          | 1            |

# Midface Plates t = 0.5 mm

## Straight Plates


Compatible with screws Ø 1.2/1.5/1.8

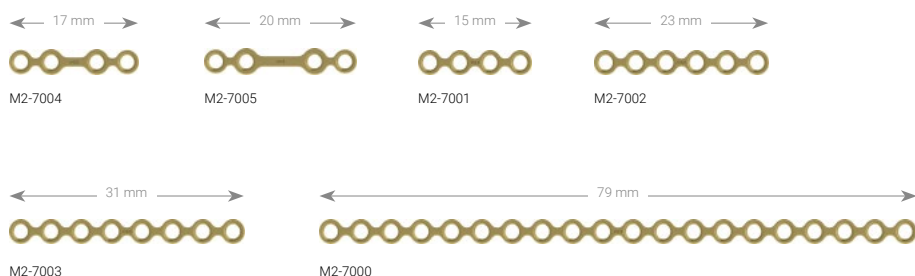Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm

| Art. No. | STERILE  | Description | Bar  | Holes | Pieces / Pkg |
|----------|----------|-------------|------|-------|--------------|
| M2-7000  | M2-7000S | straight    |      | 20    | 1            |
| M2-7001  | M2-7001S | straight    |      | 4     | 1            |
| M2-7002  | M2-7002S | straight    |      | 6     | 1            |
| M2-7003  | M2-7003S | straight    |      | 8     | 1            |
| M2-7004  | M2-7004S | straight    | 6 mm | 4     | 1            |
| M2-7005  | M2-7005S | straight    | 9 mm | 4     | 1            |

## L Plates

Compatible with screws Ø 1.2/1.5/1.8

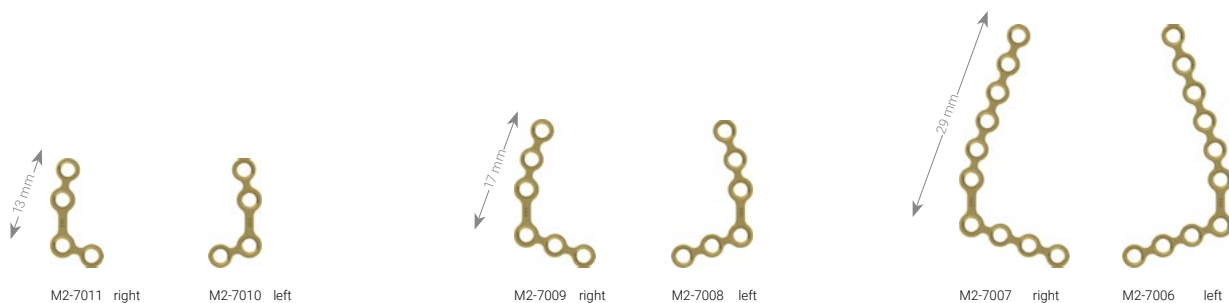
Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm

| Art. No. | STERILE  | Description  | Bar  | Holes    | Pieces / Pkg |
|----------|----------|--------------|------|----------|--------------|
| M2-7006  | M2-7006S | L left 100°  | 6 mm | 10 (4/6) | 1            |
| M2-7007  | M2-7007S | L right 100° | 6 mm | 10 (4/6) | 1            |
| M2-7008  | M2-7008S | L left 100°  | 6 mm | 6 (3/3)  | 1            |
| M2-7009  | M2-7009S | L right 100° | 6 mm | 6 (3/3)  | 1            |
| M2-7010  | M2-7010S | L left 100°  | 6 mm | 4 (2/2)  | 1            |
| M2-7011  | M2-7011S | L right 100° | 6 mm | 4 (2/2)  | 1            |

# T Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm

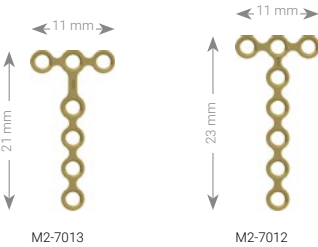


| Art. No. | STERILE  | Description | Bar  | Holes   | Pieces / Pkg |
|----------|----------|-------------|------|---------|--------------|
| M2-7012  | M2-7012S | T           |      | 8 (3/5) | 1            |
| M2-7013  | M2-7013S | T           | 6 mm | 7 (3/4) | 1            |

# Y Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm

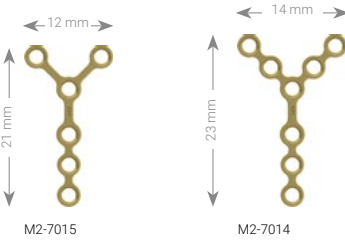


| Art. No. | STERILE  | Description | Bar  | Holes   | Pieces / Pkg |
|----------|----------|-------------|------|---------|--------------|
| M2-7014  | M2-7014S | Y           | 6 mm | 8 (5/3) | 1            |
| M2-7015  | M2-7015S | Y           | 6 mm | 6 (3/3) | 1            |

# X Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm

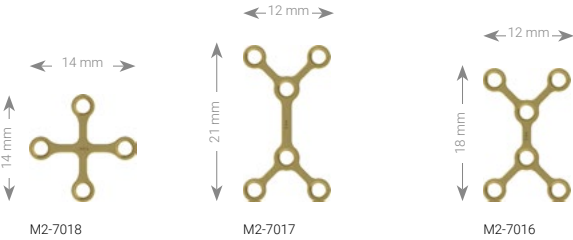


| Art. No. | STERILE  | Description | Bar  | Holes | Pieces / Pkg |
|----------|----------|-------------|------|-------|--------------|
| M2-7016  | M2-7016S | X           | 6 mm | 6     | 1            |
| M2-7017  | M2-7017S | X           | 9 mm | 6     | 1            |
| M2-7018  | M2-7018S | X           |      | 4     | 1            |

## H Plate

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm

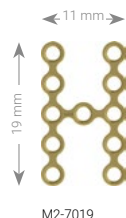

M2-7019

| Art. No. | STERILE  | Description | Bar  | Holes | Pieces / Pkg |
|----------|----------|-------------|------|-------|--------------|
| M2-7019  | M2-7019S | H           | 4 mm | 11    | 1            |

## Orbital Rim Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm



M2-7021


M2-7020

| Art. No. | STERILE  | Description | Holes | Pieces / Pkg |
|----------|----------|-------------|-------|--------------|
| M2-7020  | M2-7020S | curved      | 12    | 1            |
| M2-7021  | M2-7021S | curved      | 6     | 1            |

## Grid Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.5 mm

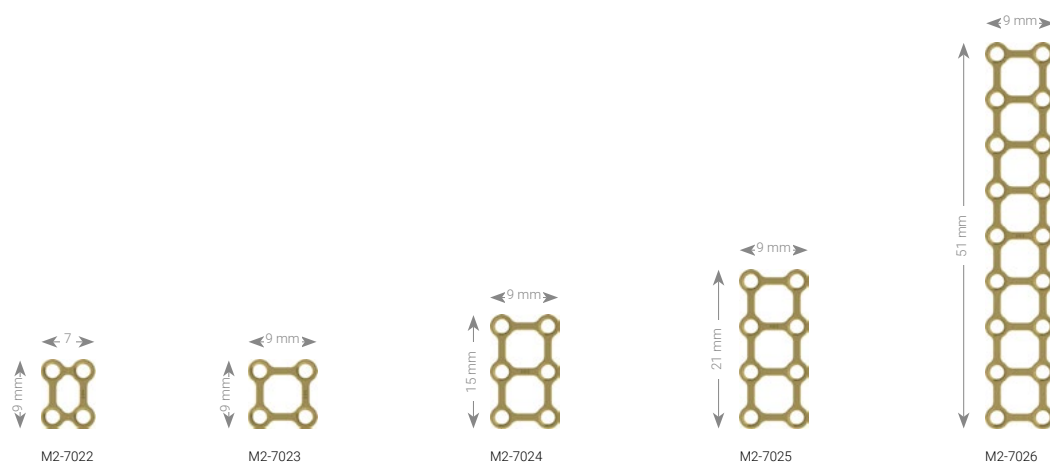

M2-7022

M2-7023

M2-7024

M2-7025


M2-7026

| Art. No. | STERILE  | Description | Holes    | Pieces / Pkg |
|----------|----------|-------------|----------|--------------|
| M2-7022  | M2-7022S | H           | 4 (2x2)  | 1            |
| M2-7023  | M2-7023S | H           | 4 (2x2)  | 1            |
| M2-7024  | M2-7024S | H           | 6 (3x2)  | 1            |
| M2-7025  | M2-7025S | H           | 8 (4x2)  | 1            |
| M2-7026  | M2-7026S | H           | 18 (9x2) | 1            |

# Midface Plates t = 0.6 / 0.7 mm

## Straight Plates

Compatible with screws Ø 1.2/1.5/1.8

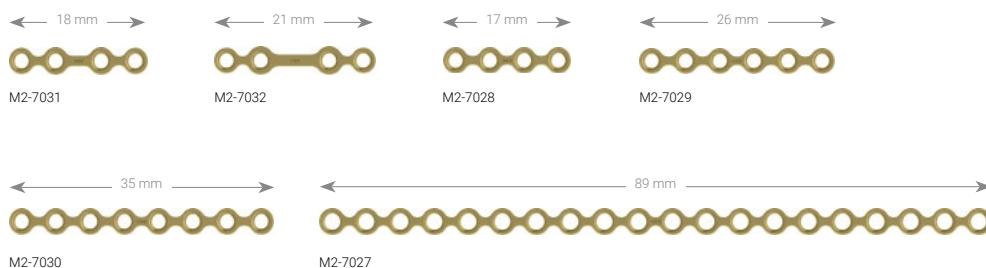Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

| Art. No. | <b>STERILE</b> | Description | Bar  | Holes | Pieces / Pkg |
|----------|----------------|-------------|------|-------|--------------|
| M2-7027  | M2-7027S       | straight    |      | 20    | 1            |
| M2-7028  | M2-7028S       | straight    |      | 4     | 1            |
| M2-7029  | M2-7029S       | straight    |      | 6     | 1            |
| M2-7030  | M2-7030S       | straight    |      | 8     | 1            |
| M2-7031  | M2-7031S       | straight    | 6 mm | 4     | 1            |
| M2-7032  | M2-7032S       | straight    | 9 mm | 4     | 1            |

## L Plates

Compatible with screws Ø 1.2/1.5/1.8


Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

| Art. No. | <b>STERILE</b> | Description  | Bar    | Holes    | Pieces / Pkg |
|----------|----------------|--------------|--------|----------|--------------|
| M2-7033  | M2-7033S       | L left 100°  | 7.5 mm | 10 (4/6) | 1            |
| M2-7034  | M2-7034S       | L right 100° | 7.5 mm | 10 (4/6) | 1            |
| M2-7035  | M2-7035S       | L left 100°  | 7.5 mm | 6 (3/3)  | 1            |
| M2-7036  | M2-7036S       | L right 100° | 7.5 mm | 6 (3/3)  | 1            |
| M2-7037  | M2-7037S       | L left 100°  | 7.5 mm | 4 (2/2)  | 1            |
| M2-7038  | M2-7038S       | L right 100° | 7.5 mm | 4 (2/2)  | 1            |

## T Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

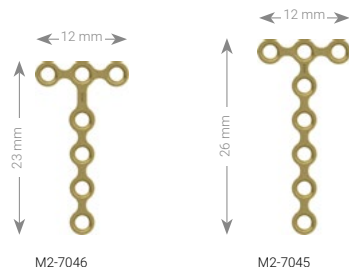


| Art. No. | STERILE  | Description | Bar  | Holes   | Pieces / Pkg |
|----------|----------|-------------|------|---------|--------------|
| M2-7045  | M2-7045S | T           |      | 8 (3/5) | 1            |
| M2-7046  | M2-7046S | T           | 6 mm | 7 (3/4) | 1            |

## Y Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

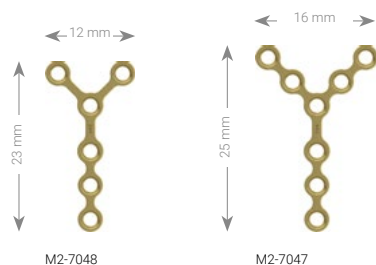


| Art. No. | STERILE  | Description | Bar  | Holes   | Pieces / Pkg |
|----------|----------|-------------|------|---------|--------------|
| M2-7047  | M2-7047S | Y           | 6 mm | 8 (5/3) | 1            |
| M2-7048  | M2-7048S | Y           | 6 mm | 6 (3/3) | 1            |

## Z Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

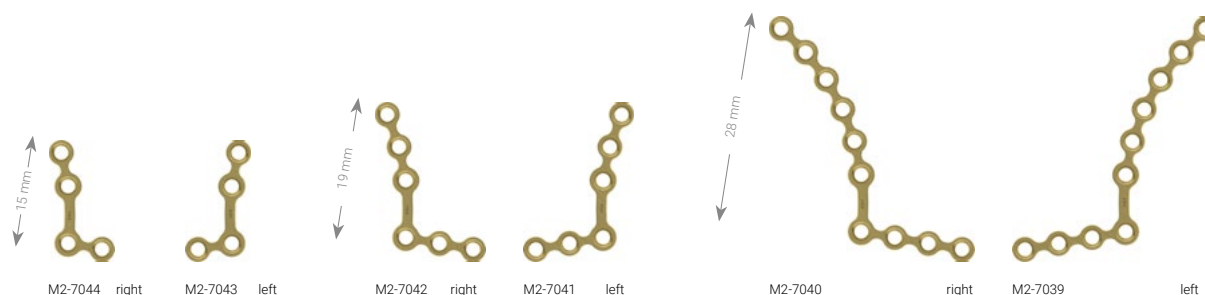


| Art. No. | STERILE  | Description  | Bar    | Holes    | Pieces / Pkg |
|----------|----------|--------------|--------|----------|--------------|
| M2-7039  | M2-7039S | Z left 100°  | 7.5 mm | 10 (4/6) | 1            |
| M2-7040  | M2-7040S | Z right 100° | 7.5 mm | 10 (4/6) | 1            |
| M2-7041  | M2-7041S | Z left 100°  | 7.5 mm | 6 (3/3)  | 1            |
| M2-7042  | M2-7042S | Z right 100° | 7.5 mm | 6 (3/3)  | 1            |
| M2-7043  | M2-7043S | Z left 100°  | 7.5 mm | 4 (2/2)  | 1            |
| M2-7044  | M2-7044S | Z right 100° | 7.5 mm | 4 (2/2)  | 1            |

## X Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

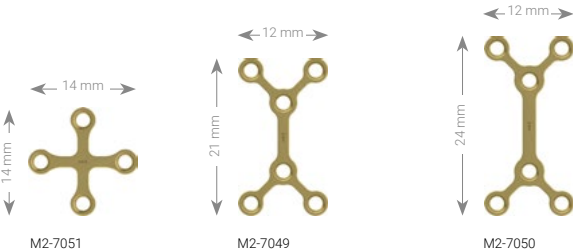


| Art. No. | STERILE  | Description | Bar   | Holes | Pieces / Pkg |
|----------|----------|-------------|-------|-------|--------------|
| M2-7049  | M2-7049S | X           | 9 mm  | 6     | 1            |
| M2-7050  | M2-7050S | X           | 12 mm | 6     | 1            |
| M2-7051  | M2-7051S | X           |       | 4     | 1            |

## H Plate

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

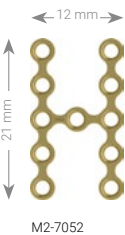


| Art. No. | STERILE  | Description | Bar    | Holes | Pieces / Pkg |
|----------|----------|-------------|--------|-------|--------------|
| M2-7052  | M2-7052S | H           | 4.5 mm | 11    | 1            |

## Orbital Rim Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.6 mm

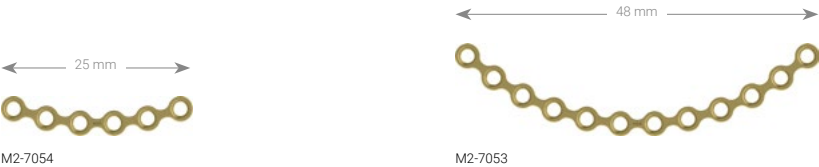


| Art. No. | STERILE  | Description | Holes | Pieces / Pkg |
|----------|----------|-------------|-------|--------------|
| M2-7053  | M2-7053S | curved      | 12    | 1            |
| M2-7054  | M2-7054S | curved      | 6     | 1            |

## Grid Plates

Compatible with screws  $\varnothing$  1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 0.7 mm

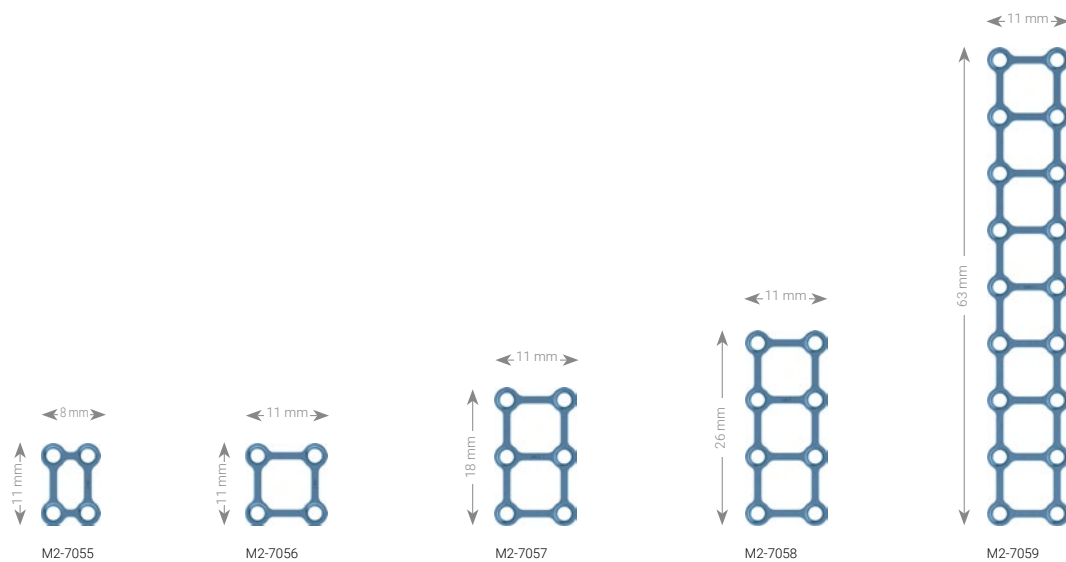


| Art. No. | STERILE  | Description | Bar          | Holes    | Pieces / Pkg |
|----------|----------|-------------|--------------|----------|--------------|
| M2-7055  | M2-7055S | H           | 4.5 x 7.5 mm | 4 (2x2)  | 1            |
| M2-7056  | M2-7056S | H           | 7.5 x 7.5 mm | 4 (2x2)  | 1            |
| M2-7057  | M2-7057S | H           | 7.5 x 7.5 mm | 6 (3x2)  | 1            |
| M2-7058  | M2-7058S | H           | 7.5 x 7.5 mm | 8 (4x2)  | 1            |
| M2-7059  | M2-7059S | H           | 7.5 x 7.5 mm | 18 (9x2) | 1            |

# Midface Plates t = 1.0 mm

## Straight Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

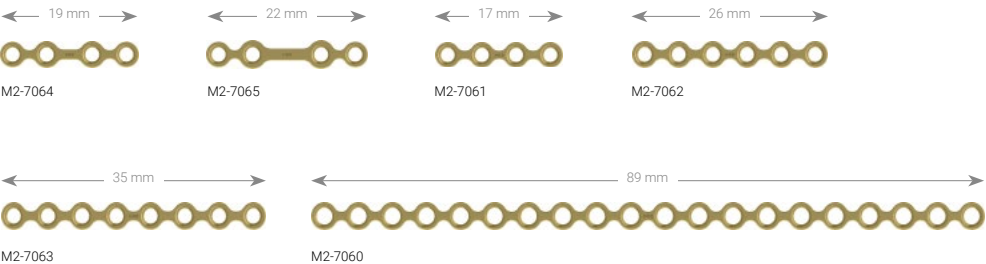


| Art. No. | STERILE  | Description | Bar  | Holes | Pieces / Pkg |
|----------|----------|-------------|------|-------|--------------|
| M2-7060  | M2-7060S | straight    |      | 20    | 1            |
| M2-7061  | M2-7061S | straight    |      | 4     | 1            |
| M2-7062  | M2-7062S | straight    |      | 6     | 1            |
| M2-7063  | M2-7063S | straight    |      | 8     | 1            |
| M2-7064  | M2-7064S | straight    | 6 mm | 4     | 1            |
| M2-7065  | M2-7065S | straight    | 9 mm | 4     | 1            |

## L Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

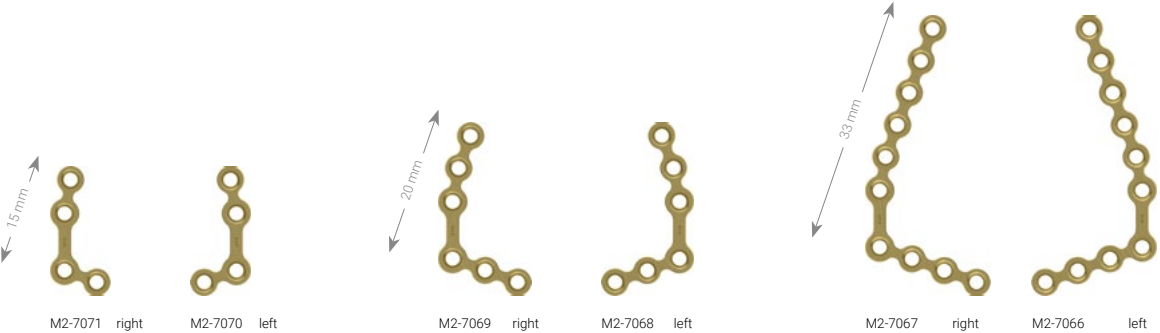


| Art. No. | STERILE  | Description  | Bar    | Holes    | Pieces / Pkg |
|----------|----------|--------------|--------|----------|--------------|
| M2-7066  | M2-7066S | L left 100°  | 7.5 mm | 10 (4/6) | 1            |
| M2-7067  | M2-7067S | L right 100° | 7.5 mm | 10 (4/6) | 1            |
| M2-7068  | M2-7068S | L left 100°  | 7.5 mm | 6 (3/3)  | 1            |
| M2-7069  | M2-7069S | L right 100° | 7.5 mm | 6 (3/3)  | 1            |
| M2-7070  | M2-7070S | L left 100°  | 7.5 mm | 4 (2/2)  | 1            |
| M2-7071  | M2-7071S | L right 100° | 7.5 mm | 4 (2/2)  | 1            |

## T Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

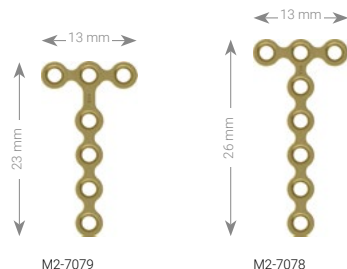


| Art. No. | STERILE  | Description | Bar  | Holes   | Pieces / Pkg |
|----------|----------|-------------|------|---------|--------------|
| M2-7078  | M2-7078S | T           |      | 8 (3/5) | 1            |
| M2-7079  | M2-7079S | T           | 6 mm | 7 (3/4) | 1            |

## Y Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

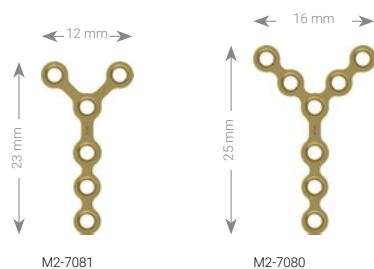


| Art. No. | STERILE  | Description | Bar  | Holes   | Pieces / Pkg |
|----------|----------|-------------|------|---------|--------------|
| M2-7080  | M2-7080S | Y           | 6 mm | 8 (5/3) | 1            |
| M2-7081  | M2-7081S | Y           | 6 mm | 6 (3/3) | 1            |

## Z Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

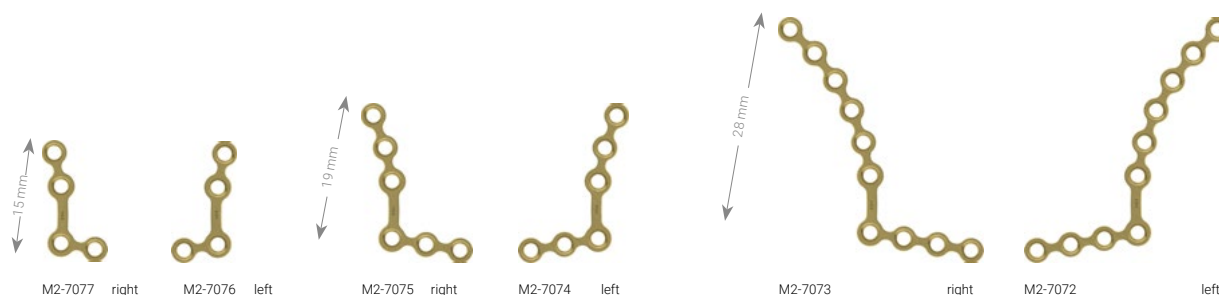


| Art. No. | STERILE  | Description  | Bar    | Holes    | Pieces / Pkg |
|----------|----------|--------------|--------|----------|--------------|
| M2-7072  | M2-7072S | Z left 100°  | 7.5 mm | 10 (4/6) | 1            |
| M2-7073  | M2-7073S | Z right 100° | 7.5 mm | 10 (4/6) | 1            |
| M2-7074  | M2-7074S | Z left 100°  | 7.5 mm | 6 (3/3)  | 1            |
| M2-7075  | M2-7075S | Z right 100° | 7.5 mm | 6 (3/3)  | 1            |
| M2-7076  | M2-7076S | Z left 100°  | 7.5 mm | 4 (2/2)  | 1            |
| M2-7077  | M2-7077S | Z right 100° | 7.5 mm | 4 (2/2)  | 1            |

## X Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

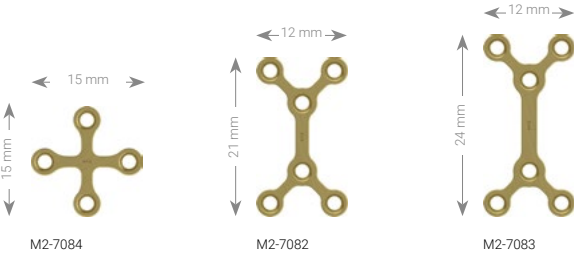


| Art. No. | STERILE  | Description | Bar   | Holes | Pieces / Pkg |
|----------|----------|-------------|-------|-------|--------------|
| M2-7082  | M2-7082S | X           | 9 mm  | 6     | 1            |
| M2-7083  | M2-7083S | X           | 12 mm | 6     | 1            |
| M2-7084  | M2-7084S | X           |       | 4     | 1            |

## H Plate

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

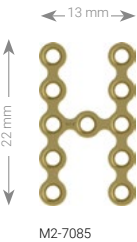


| Art. No. | STERILE  | Description | Bar    | Holes | Pieces / Pkg |
|----------|----------|-------------|--------|-------|--------------|
| M2-7085  | M2-7085S | H           | 4.5 mm | 11    | 1            |

## Orbital Rim Plates

Compatible with screws Ø 1.2/1.5/1.8

Material: cpTi (ASTM F67)  
Plate thickness: 1.0 mm

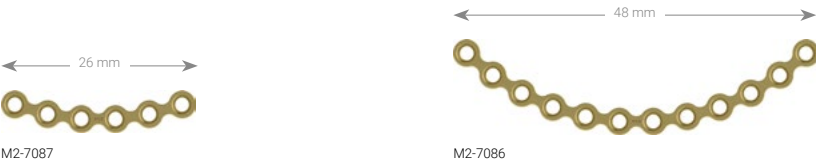

| Art. No. | STERILE  | Description | Holes | Pieces / Pkg |
|----------|----------|-------------|-------|--------------|
| M2-7086  | M2-7086S | curved      |       | 121          |
| M2-7087  | M2-7087S | curved      |       | 61           |

## 1.2 Cortical Screws, HexaDrive 4

Material: Ti6Al4V (ASTM F136)



| Length | Art. No.     | STERILE       | Pieces / Pkg |
|--------|--------------|---------------|--------------|
| 3 mm   | M2-5214.03/1 | M2-5214.03/1S | 1            |
| 4 mm   | M2-5214.04/1 | M2-5214.04/1S | 1            |
| 5 mm   | M2-5214.05/1 | M2-5214.05/1S | 1            |
| 6 mm   | M2-5214.06/1 | M2-5214.06/1S | 1            |
| 7 mm   | M2-5214.07/1 | M2-5214.07/1S | 1            |
| 8 mm   | M2-5214.08/1 | M2-5214.08/1S | 1            |
| 9 mm   | M2-5214.09/1 | M2-5214.09/1S | 1            |
| 11 mm  | M2-5214.11/1 | M2-5214.11/1S | 1            |

| Art. No.   | Pieces / Pkg |
|------------|--------------|
| M2-5214.03 | 5            |
| M2-5214.04 | 5            |
| M2-5214.05 | 5            |
| M2-5214.06 | 5            |
| M2-5214.07 | 5            |
| M2-5214.08 | 5            |
| M2-5214.09 | 5            |
| M2-5214.11 | 5            |

## 1.5 Cortical Screws, HexaDrive 4

Material: Ti6Al4V (ASTM F136)



| Length | Art. No.     | STERILE       | Pieces / Pkg |
|--------|--------------|---------------|--------------|
| 3 mm   | M2-5224.03/1 | M2-5224.03/1S | 1            |
| 4 mm   | M2-5224.04/1 | M2-5224.04/1S | 1            |
| 5 mm   | M2-5224.05/1 | M2-5224.05/1S | 1            |
| 6 mm   | M2-5224.06/1 | M2-5224.06/1S | 1            |
| 7 mm   | M2-5224.07/1 | M2-5224.07/1S | 1            |
| 8 mm   | M2-5224.08/1 | M2-5224.08/1S | 1            |
| 9 mm   | M2-5224.09/1 | M2-5224.09/1S | 1            |
| 11 mm  | M2-5224.11/1 | M2-5224.11/1S | 1            |

| Art. No.   | Pieces / Pkg |
|------------|--------------|
| M2-5224.03 | 5            |
| M2-5224.04 | 5            |
| M2-5224.05 | 5            |
| M2-5224.06 | 5            |
| M2-5224.07 | 5            |
| M2-5224.08 | 5            |
| M2-5224.09 | 5            |
| M2-5224.11 | 5            |

## 1.5 SpeedTip Screws, HexaDrive 4

Material: Ti6Al4V (ASTM F136)



| Length | Art. No.     | STERILE       | Pieces / Pkg |
|--------|--------------|---------------|--------------|
| 4 mm   | M2-5223.04/1 | M2-5223.04/1S | 1            |
| 5 mm   | M2-5223.05/1 | M2-5223.05/1S | 1            |
| 6 mm   | M2-5223.06/1 | M2-5223.06/1S | 1            |
| 7 mm   | M2-5223.07/1 | M2-5223.07/1S | 1            |
| 8 mm   | M2-5223.08/1 | M2-5223.08/1S | 1            |
| 9 mm   | M2-5223.09/1 | M2-5223.09/1S | 1            |

| Art. No.   | Pieces / Pkg |
|------------|--------------|
| M2-5223.04 | 5            |
| M2-5223.05 | 5            |
| M2-5223.06 | 5            |
| M2-5223.07 | 5            |
| M2-5223.08 | 5            |
| M2-5223.09 | 5            |

## 1.8 Cortical Screws, HexaDrive 4

Material: Ti6Al4V (ASTM F136)



| Length | Art. No.     | STERILE       | Pieces / Pkg |
|--------|--------------|---------------|--------------|
| 3 mm   | M2-5234.03/1 | M2-5234.03/1S | 1            |
| 4 mm   | M2-5234.04/1 | M2-5234.04/1S | 1            |
| 5 mm   | M2-5234.05/1 | M2-5234.05/1S | 1            |
| 6 mm   | M2-5234.06/1 | M2-5234.06/1S | 1            |
| 7 mm   | M2-5234.07/1 | M2-5234.07/1S | 1            |
| 8 mm   | M2-5234.08/1 | M2-5234.08/1S | 1            |
| 9 mm   | M2-5234.09/1 | M2-5234.09/1S | 1            |
| 11 mm  | M2-5234.11/1 | M2-5234.11/1S | 1            |

| Art. No.   | Pieces / Pkg |
|------------|--------------|
| M2-5234.03 | 5            |
| M2-5234.04 | 5            |
| M2-5234.05 | 5            |
| M2-5234.06 | 5            |
| M2-5234.07 | 5            |
| M2-5234.08 | 5            |
| M2-5234.09 | 5            |
| M2-5234.11 | 5            |

Twist Drills Ø 1.0 mm (Core Hole 1.2 Screws)



M2-3012



M2-3022



M2-3032



M2-3042



M2-3052



M2-3062



M2-3382



M2-3392

| Art. No. | <b>STERILE</b> | Description             | Stop  | Length | Shaft End       | Pieces / Pkg |
|----------|----------------|-------------------------|-------|--------|-----------------|--------------|
| M2-3012  | M2-3012S       |                         | 5 mm  | 35 mm  | Dental          | 1            |
| M2-3022  | M2-3022S       |                         | 5 mm  | 48 mm  | Stryker J-Latch | 1            |
| M2-3032  | M2-3032S       |                         | 7 mm  | 37 mm  | Dental          | 1            |
| M2-3042  | M2-3042S       |                         | 7 mm  | 50 mm  | Stryker J-Latch | 1            |
| M2-3052  | M2-3052S       |                         | 25 mm | 55 mm  | Dental          | 1            |
| M2-3062  | M2-3062S       |                         | 25 mm | 68 mm  | Stryker J-Latch | 1            |
| M2-3382  | M2-3382S       | for drill guide M2-2202 | 25 mm | 99 mm  | Dental          | 1            |
| M2-3392  | M2-3392S       | for drill guide M2-2202 | 25 mm | 112 mm | Stryker J-Latch | 1            |

Twist Drills Ø 1.2 mm (Gliding Hole 1.2 Screws)



M2-3072



M2-3322



M2-3082



M2-3332

| Art. No. | <b>STERILE</b> | Description             | Stop  | Length | Shaft End       | Pieces / Pkg |
|----------|----------------|-------------------------|-------|--------|-----------------|--------------|
| M2-3072  | M2-3072S       |                         | 25 mm | 55 mm  | Dental          | 1            |
| M2-3082  | M2-3082S       |                         | 25 mm | 68 mm  | Stryker J-Latch | 1            |
| M2-3322  | M2-3322S       | for drill guide M2-2202 | 25 mm | 99 mm  | Dental          | 1            |
| M2-3332  | M2-3332S       | for drill guide M2-2202 | 25 mm | 112 mm | Stryker J-Latch | 1            |

## Twist Drills Ø 1.2 mm (Core Hole 1.5 Screws)



M2-3122



M2-3132



M2-3142



M2-3152



M2-3162



M2-3172



M2-3402



M2-3412

| Art. No. | STERILE  | Description             | Stop  | Length | Shaft End       | Pieces / Pkg |
|----------|----------|-------------------------|-------|--------|-----------------|--------------|
| M2-3122  | M2-3122S |                         | 5 mm  | 35 mm  | Dental          | 1            |
| M2-3132  | M2-3152S |                         | 5 mm  | 48 mm  | Stryker J-Latch | 1            |
| M2-3142  | M2-3132S |                         | 7 mm  | 37 mm  | Dental          | 1            |
| M2-3152  | M2-3162S |                         | 7 mm  | 50 mm  | Stryker J-Latch | 1            |
| M2-3162  | M2-3142S |                         | 25 mm | 55 mm  | Dental          | 1            |
| M2-3172  | M2-3172S |                         | 25 mm | 68 mm  | Stryker J-Latch | 1            |
| M2-3402  | M2-3402S | for drill guide M2-2202 | 25 mm | 99 mm  | Dental          | 1            |
| M2-3412  | M2-3412S | for drill guide M2-2202 | 25 mm | 112 mm | Stryker J-Latch | 1            |

## Twist Drills Ø 1.5 mm (Gliding Hole 1.5 Screws)



M2-3182



M2-3342



M2-3192



M2-3352

| Art. No. | STERILE  | Description             | Stop  | Length | Shaft End       | Pieces / Pkg |
|----------|----------|-------------------------|-------|--------|-----------------|--------------|
| M2-3182  | M2-3182S |                         | 25 mm | 55 mm  | Dental          | 1            |
| M2-3192  | M2-3192S |                         | 25 mm | 68 mm  | Stryker J-Latch | 1            |
| M2-3342  | M2-3342S | for drill guide M2-2202 | 25 mm | 99 mm  | Dental          | 1            |
| M2-3352  | M2-3352S | for drill guide M2-2202 | 25 mm | 112 mm | Stryker J-Latch | 1            |

# Twist Drills Ø 1.5 mm (Core Hole 1.8 Screws)



M2-3212



M2-3222



M2-3232



M2-3242



M2-3252



M2-3262



M2-3422



M2-3452

| Art. No. | <b>STERILE</b> | Description             | Stop  | Length | Shaft End       | Pieces / Pkg |
|----------|----------------|-------------------------|-------|--------|-----------------|--------------|
| M2-3212  | M2-3212S       |                         | 5 mm  | 35 mm  | Dental          | 1            |
| M2-3222  | M2-3222S       |                         | 5 mm  | 48 mm  | Stryker J-Latch | 1            |
| M2-3232  | M2-3232S       |                         | 7 mm  | 37 mm  | Dental          | 1            |
| M2-3242  | M2-3242S       |                         | 7 mm  | 50 mm  | Stryker J-Latch | 1            |
| M2-3252  | M2-3252S       |                         | 25 mm | 55 mm  | Dental          | 1            |
| M2-3262  | M2-3262S       |                         | 25 mm | 68 mm  | Stryker J-Latch | 1            |
| M2-3422  | M2-3422S       | for drill guide M2-2202 | 25 mm | 99 mm  | Dental          | 1            |
| M2-3452  | M2-3452S       | for drill guide M2-2202 | 25 mm | 112 mm | Stryker J-Latch | 1            |

# Twist Drills Ø 1.8 mm (Gliding Hole 1.8 Screws)



M2-3272



M2-3362



M2-3282



M2-3372

| Art. No. | <b>STERILE</b> | Description             | Stop  | Length | Shaft End       | Pieces / Pkg |
|----------|----------------|-------------------------|-------|--------|-----------------|--------------|
| M2-3272  | M2-3272S       |                         | 25 mm | 55 mm  | Dental          | 1            |
| M2-3282  | M2-3282S       |                         | 25 mm | 68 mm  | Stryker J-Latch | 1            |
| M2-3362  | M2-3362S       | for drill guide M2-2202 | 25 mm | 99 mm  | Dental          | 1            |
| M2-3372  | M2-3372S       | for drill guide M2-2202 | 25 mm | 112 mm | Stryker J-Latch | 1            |

## Drill Guide



| Art. No. | System Size | Length | Pieces / Pkg |
|----------|-------------|--------|--------------|
| M2-2202  | 1.2–1.8     | 164 mm | 1            |

## Depth Gauge



M2-2250

| Art. No. | System Size | Length | Pieces / Pkg |
|----------|-------------|--------|--------------|
| M2-2250  | 1.2–2.3     | 153 mm | 1            |

## Screwdriver Handles



M2-2001



M2-2003

| Art. No. | Description                             | Length | Pieces / Pkg |
|----------|-----------------------------------------|--------|--------------|
| M2-2001  | type 2 (hand-driven small, AO coupling) | 121 mm | 1            |
| M2-2003  | type 1 (finger-driven, AO coupling)     | 121 mm | 1            |

## Screwdriver Blade



1:1

| Art. No. | Interface | Description  | Length | Pieces / Pkg |
|----------|-----------|--------------|--------|--------------|
| M2-2004  | HD4       | self-holding | 80 mm  | 1            |

Plate and Screw Holding Forceps




M-2009




M-2019

| Art. No. | Description   | Length | Pieces / Pkg |
|----------|---------------|--------|--------------|
| M-2009   | angled, small | 150 mm | 1            |
| M-2019   | angled, large | 200 mm | 1            |

Plate Cutting Pliers

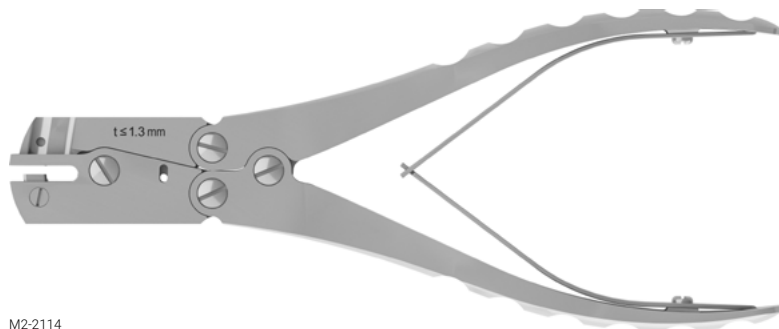

M2-2114

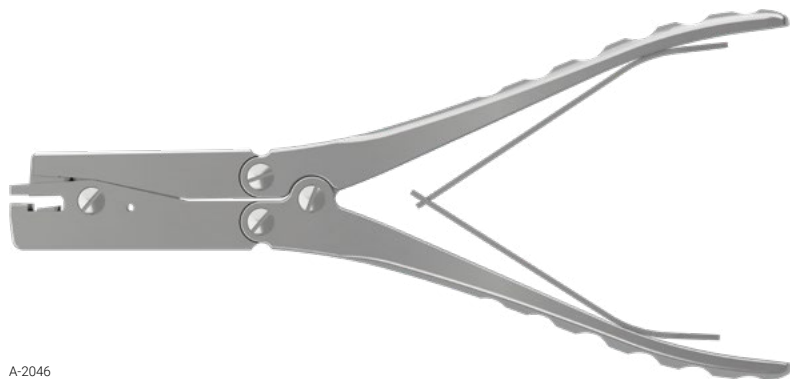


A-2046

| Art. No. | System Size | Description | Length | Pieces / Pkg |
|----------|-------------|-------------|--------|--------------|
| M2-2114  |             | t ≤ 1.3 mm  | 204 mm | 1            |
| A-2046   | 1.2–2.8     |             | 207 mm | 1            |

## Mesh Cutting Pliers



M2-2870

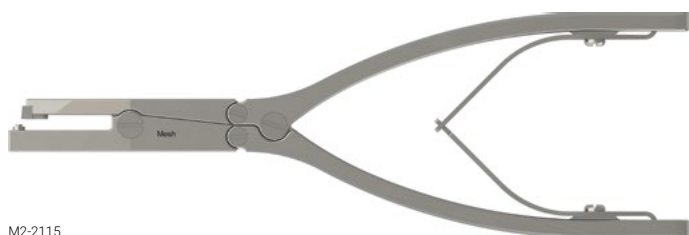

M2-2115

| Art. No. | Description | Length | Pieces / Pkg |
|----------|-------------|--------|--------------|
| M2-2115  |             | 180 mm | 1            |
| M2-2870  | curved left | 127 mm | 1            |

## Plate Bending Pliers



M2-2000



M2-2002



M2-2012

| Art. No. | System Size | Description | Length | Pieces / Pkg |
|----------|-------------|-------------|--------|--------------|
| M2-2000  |             | flat        | 152 mm | 1            |
| M2-2002  | 1.2-1.8     |             | 152 mm | 1            |
| M2-2012  | 1.2-1.8     | with pin    | 140 mm | 1            |

Orbital Retractors



M2-2121



M2-2122



M2-2123


| Art. No. | Description | Length | Pieces / Pkg |
|----------|-------------|--------|--------------|
| M2-2121  | right       | 200 mm | 1            |
| M2-2122  | left        | 200 mm | 1            |
| M2-2123  |             | 200 mm | 1            |

Elevator and Mesh Forming Instrument



| Art. No. | Length | Pieces / Pkg |
|----------|--------|--------------|
| M2-2872  | 190 mm | 1            |

Midface Model

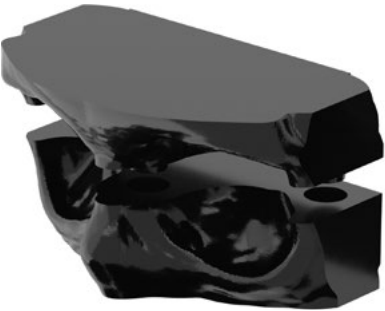

| Art. No. | Description | Length | Pieces / Pkg |
|----------|-------------|--------|--------------|
| M2-2874  | two pieces  | 101 mm | 1            |



R\_MIDFACE2-02000001\_v0 / 2025-05, Medartis AG, Switzerland. All technical data subject to alteration.

## MANUFACTURER & HEADQUARTERS

Medartis AG | Hochbergerstrasse 60E | 4057 Basel / Switzerland

P +41 61 633 34 34 | F +41 61 633 34 00 | [www.medartis.com](http://www.medartis.com)

## SUBSIDIARIES

Australia | Austria | Brazil | France | Germany | Japan | Mexico | New Zealand | Poland | Spain | UK | USA

For detailed information regarding our subsidiaries and distributors, please visit [www.medartis.com](http://www.medartis.com)



Disclaimer: This information is intended to demonstrate the Medartis portfolio of medical devices. A surgeon must always rely on her or his own professional clinical judgement when deciding whether to use a particular product when treating a particular patient. Medartis is not giving any medical advice. The devices may not be available in all countries due to registration and/or medical practices. For further questions, please contact your Medartis representative ([www.medartis.com](http://www.medartis.com)). This information contains products with CE and/or UKCA marking. All pictures shown are for illustration purposes only and may not be an exact representation of the product. For US only: Federal law restricts this device to sale by or on the order of a physician.

© Medartis 2025. Everything herein is protected by copyright, trademarks and other intellectual property rights, as applicable, owned by or licensed to Medartis or its affiliates unless otherwise indicated. It is forbidden to redistribute, duplicate or disclose anything herein, in whole or in part, without the prior written consent of Medartis.